

A-Trac **Dynamic Brace**™ Patent Pending

The Anti-Translational Cruciate Brace
by WoundWear, Inc.



Premier Edition


- ◆ Improve your treatment outcomes
- ◆ Provide an effective treatment alternative for dogs not undergoing surgery
- ◆ Enhance client satisfaction

Treatment Issues:

- There is a 30 to 60% incidence of contra lateral cruciate injury.
- Frequent and unnecessary re-injuries following cruciate surgical repair as a result of a stifle joint which is left unprotected post-operatively.
- Lack of effective treatment alternatives for poor surgical candidates or clients who can't afford surgery.
- Delayed treatment due to wait time for surgery, or "off loading" of the affected limb following cruciate injury or surgery, causes atrophy of thigh muscles and delays the reestablishment of proprioception.
- Arthritis progresses as long as stifle joint instability is allowed to continue.

Brace Function & Features:

- The brace's tri-laminated material provides comfort, elasticity, and conformability to the dog's leg enhancing support and proprioception.
- Sagittal plane stability is provided by patent pending tension rods tethered at an appropriate distance above and below the stifle joint on each side of the brace. Two supplied sets of tension rods provide for different ranges of motion against variable resistance. This imparts continuous rehabilitation & allows the veterinarian individualized control over his rehabilitation protocol.
- Anterior translation of the tibia on the femur is prevented by a buttress or blocking mechanism between the two tension rods. Stabilization of the stifle diminishes pain, arthritis progression, and allows for normal weight bearing on the affected limb. The dog does not perceive a need to off-load weight to the unaffected limb.

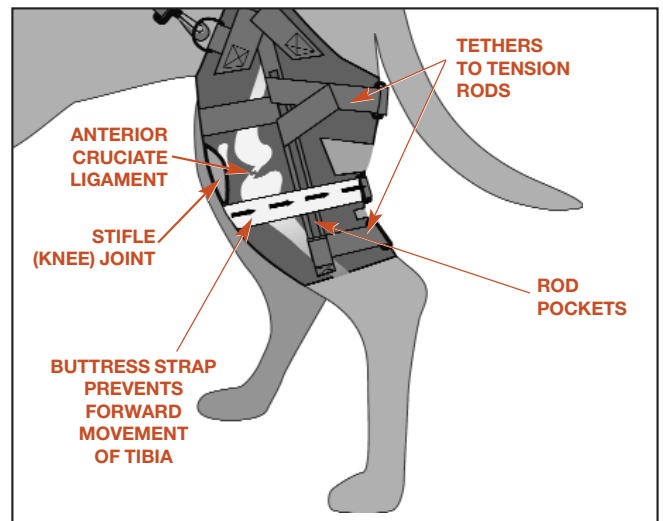

DON'T LEAVE THE SUCCESS OF YOUR TREATMENT TO CHANCE.

Prescribe the A-TraC Dynamic Brace™:

- Immediately after a cruciate injury has been diagnosed
- Continue brace treatment after surgery, or as part of a non-surgical treatment regimen



To order by phone, or request product literature from the United States, its principalities or Canada, please use our toll free number: **800-443-4055**

Outside the United States, please call: **970-679-8636**

For Internet orders, or product information, please see our web site at www.woundwear.com

For assistance with our products, call: **847-634-1700** WoundWear, Inc. P.O. Box 6008 Buffalo Grove, Illinois 60089